

e-Ordering

You can efficiently order LifeVest, or renew an existing medical order, through the LifeVest Network's e-Ordering interface. Once submitted, the LifeVest medical order is transmitted directly to ZOLL for processing.

The e-Ordering functionality allows for coordination of efforts of the clinical staff (e.g., upon direction from a physician, a staff member may enter information into the online medical order form, attach proper documentation, and forward to the prescriber for review and electronic signature.)

The screenshot shows the ZOLL LifeVest e-Ordering interface. At the top, there is a search bar for the patient's last name and a login section. The main form contains several sections: Prescribing Physician (with a dropdown for NP), Patient Name (First Name, MI, Last Name), Date of Birth (Year, Month, Day), Primary Insurance (Primary and Secondary), Reason for LifeVest (with a dropdown for Set Primary Reason), and Upload Supporting Documentation (with a 'Choose File' button and an 'Attach File' button). Below these are fields for Estimated Length of Need and Estimated Start Date. There are also sections for VF Rate Threshold (BPM and BPM) and Emergency Injuries (Shock 1 through Shock 5). A 'Reset to Defaults' button is present. A comments field is also available. At the bottom, there are two main options: 'e-Order LifeVest' and 'Fax Order LifeVest'. The 'e-Order LifeVest' option lists: 'Electronically submit medical order', 'ZOLL will begin processing', and 'Order will be simultaneously sent electronically to prescriber for e-signature'. The 'Fax Order LifeVest' option lists: 'Generates PDF of medical order with fax cover page' and 'PDF must be printed, signed by prescriber, and faxed to ZOLL'. An 'OR' is placed between these two options. A 'Save and Resume Later' button is at the bottom center. A small note at the bottom left says '* Denotes Required Field'.

Contact **ZOLL** 24 hours a day, 7 days a week



For LifeVest customer support, technical support or medical orders, please call **800.543.3267**



LifeVest medical orders and supporting documentation can be faxed to **866.567.7615**



For more LifeVest information, please visit www.zoll.com or email LifeVest.Info@zoll.com



Access LifeVest patient data and e-Order through LifeVestNetwork.Zoll.com

*When the rhythm meets heart rate and morphology criteria, or through a patient initiated recording.



Proven Therapy for
Patients at Risk of SCD
+ so much more™



The LifeVest Network

Patient Data Management System

The LifeVest® wearable cardioverter defibrillator (WCD) is a proven therapy for patients at risk of Sudden Cardiac Death (SCD) and provides a connection to your patient's changing cardiac condition through the LifeVest Network.

The LifeVest Network is ZOLL's online patient data management system that allows clinicians to monitor patient data downloaded from a patient's LifeVest WCD. The LifeVest Network enables you to have access to the most current information about your patient during their recovery after a cardiac event.


Data captured by the patient's LifeVest device is transmitted to the secure LifeVest Network website where it is archived for review. Information on the LifeVest Network can be accessed anytime from almost anywhere.

ECG Review

ECG recordings automatically recorded by LifeVest or manually captured by the patient provide additional information that can be used to decide future care paths.

ECG information captured by LifeVest can help you diagnose:


- Sustained VT/VF
- Non-sustained VT
- Atrial arrhythmias/SVT*
- Severe bradycardia/Asystole

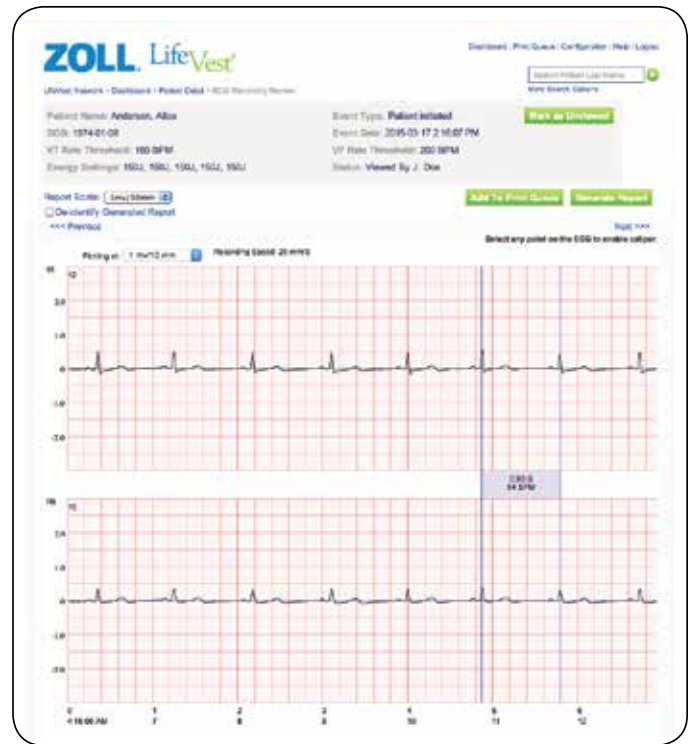

Trends

LifeVest Network Trends provides you with continuous insight into your patients' health status and response to therapy during the early period after a cardiac event, offering information that you can use to optimize their care.

Information about your patient is trended over time, indicating when a change in their status may require your action.

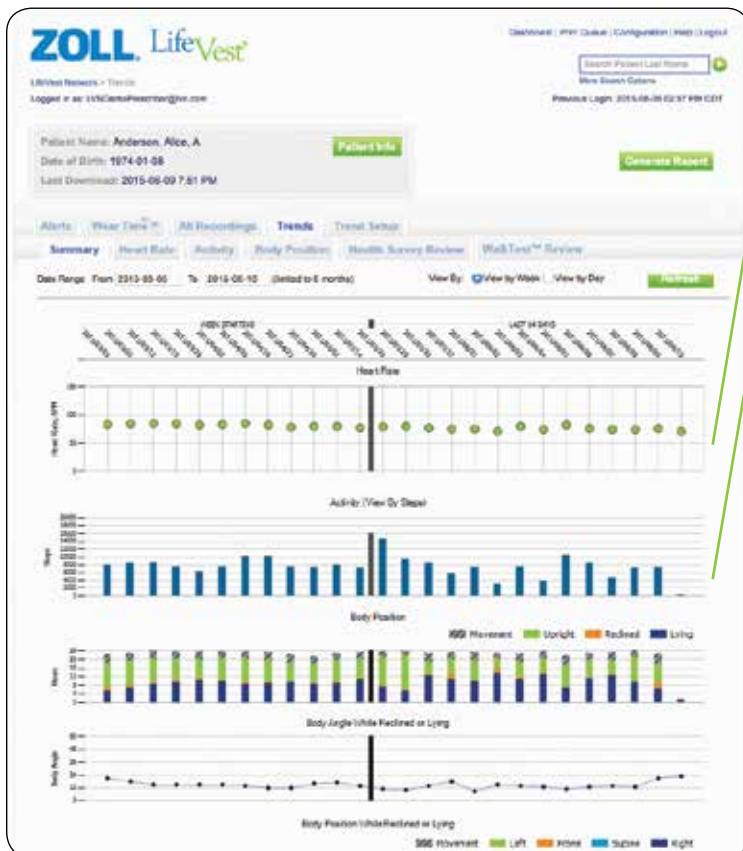

Heart Rate

Heart rate displayed as weekly, daily, or five-minute averages

Activity Measures

- Patient activity (steps taken)
- Body position
- Body angle while at rest

Health Survey

The optional Health Survey allows you to select up to 12 questions that your patients can answer through the LifeVest's touchscreen display.

WalkTest™

You can prescribe an optional six minute walk test. The LifeVest monitor will guide your patient through the test using the LifeVest's touchscreen display. Answers to pre/post walk questions, distance walked, and heart rate are displayed.